Excel Xena EVO 3+4

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User manual

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Table of contents

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Table des matières

Nederlands.....4

English.....30

Deutsch.....56

Français.....82

Nederlands

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English

Table of contents

General introduction	32
1. Identification	34
2. General information	35
2.1 Components and structure of your mobility scooter	35
2.2 General explanation of your mobility scooter	35
3. Safety instructions	36
3.1 General safety instructions	36
3.2 Warnings for safe use	36
4. Specifications	37
5. Components of your Excel Xena EVO mobility scooter	38
5.1 The dashboard	38
5.2 The tiller	39
5.3 The armrests	39
5.4 The seat	40
5.5 The freewheel	41
5.6 The battery	41
5.6.1 Charging the battery	41
5.6.2 Removing the battery	42
5.6.3 The fuse	43
5.7 The anti-tip wheels	43
5.8 The wheels	43
5.8.1 Front wheels	43
5.8.2 Rear wheels	43
6. User instructions	44
6.1 Balance	44
6.1.1 Getting in and out	44

6.2	Bending and/or reaching.....	45
6.2.1	Reaching forward.....	45
6.2.2	Reaching sideways.....	45
6.2.3	Reaching behind.....	45
6.3	Navigating kerbs and slopes.....	45
6.3.1	Going up a kerb or slope.....	45
6.3.2	Going down a kerb or slope.....	45
6.3.3	Moving surfaces.....	46
6.4	Steering.....	46
6.4.1	Steering in different situations.....	46
6.5	EMI / RFI.....	47
6.5.1	EMI / RFI recommendations.....	47
7	Transport and handling of the mobility scooter.....	48
7.1	Dismantling the scooter for transport.....	48
7.2	Reassembling the scooter for use.....	49
8	Maintenance.....	50
8.1	Tyres.....	50
8.2	Bodywork.....	50
8.3	Upholstery.....	50
8.4	Cleaning your mobility scooter.....	51
8.5	Periodic checks.....	51
9.	Troubleshooting.....	52
10.	Warranty.....	53
10.1	Warranty coverage.....	53
10.2	Warranty terms and conditions.....	53
10.2.1	Warranty terms for LiFePO4 batteries.....	55

General introduction

Before using your new Excel Xena EVO mobility scooter, both you and any potential caregiver must read and understand this user manual in full. By reading and understanding this manual, you also agree to the terms and conditions set out within it.

You are now becoming familiar with a quality product offered by Van Os Medical B.V. We thank you for placing your trust in the Excel Mobility® brand. The Excel Xena EVO mobility scooter is a true quality product. Van Os Medical B.V. is committed to continuously developing the quality and reliability of all products in its range. We therefore reserve the right to make changes to this user manual without prior notice.

It is important that you read this manual carefully before using your mobility scooter. This manual contains essential information about the safe use and maintenance of your scooter. We recommend that you keep this user manual in a safe place for future reference, ensuring safe and proper use.

The safety instructions in this manual are general guidelines and should be seen as broad recommendations. You may develop your own methods for common tasks. However, we strongly advise consulting a qualified professional to help you develop safe and effective techniques suited to your daily activities and physical capabilities.

Your new mobility scooter requires regular, scheduled maintenance. While much of this can be done by the user, we recommend having your scooter inspected at least once a year by a professional. Failure to have the scooter maintained in a timely manner will void the warranty.

Note!

In this user manual you will also find tips and warnings. These are clearly marked using symbols and styled text, as shown below.

Tip

The tips in this manual are intended to help you make the most of your Excel Xena EVO mobility scooter.

Warning

The warnings included in this manual must always be followed to prevent damage to the Excel Xena EVO mobility scooter or injury to yourself.

Fill in the details of your authorised dealer and product below:

Company name:

Address:

Phone number:

Email address:

Website:

Purchase date:


If applicable, stamp or sticker from your authorised dealer:

If there was an additional detachable frame label included with your product alongside the one on the product, you can place it below for easier reference in the future:

1. Identification

Your mobility scooter is equipped with a frame label, located on the seat column (underneath the seat) of your scooter. The frame label provides key information about your scooter, such as the product name, manufacturer name, name of the European authorised representative, Unique Device Identifier (UDI), and technical details. Below is an example of a frame label for illustration purposes. The label on your own scooter may differ. The data on your label correspond to your specific product.

Van Os Medical B.V. is the importer of this Excel Mobility product. The product is manufactured by Excel Mobility Manufacturing Co. Ltd. and is represented in Europe by SUNGO Europe B.V.

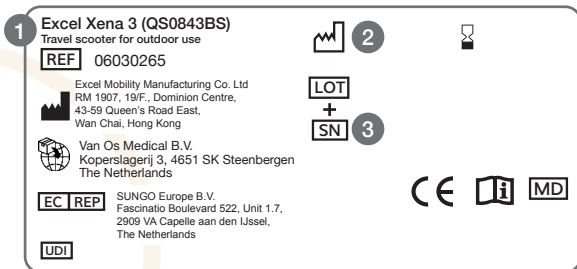1. Model name + description

The model name of your scooter starts with the brand name Excel. The brand name Excel combined with a descriptive model name forms the full model name of your scooter. When requesting information or support, it is important to have the model name, in this case Excel Xena EVO, ready. You will also find the intended use of your scooter here.

2. Date of manufacture

The date on which your mobility scooter was manufactured.

3. Serial number

Each mobility scooter has its own unique serial number. You will need this number when making technical enquiries or ordering parts for your scooter.

Warning

Your Excel mobility scooter must only be used for its intended purpose.

2. General information

2.1 Components and structure of your mobility scooter

1. Dashboard
2. Tiller
3. Basket
4. Lighting
5. Front wheel
6. Battery pack
7. Rear wheel
8. Anti-tip wheel
9. Freewheel lever
10. Seat
11. Armrest

Photo 1

Your mobility scooter is equipped with various components and parts. Familiarise yourself with these before continuing to study this manual. Designs and specifications are subject to change without prior notice.

2.2 General explanation of your mobility scooter

Your new Excel Xena EVO is a travel scooter, ideal for both indoor and outdoor use. It features full suspension (both front and rear), ensuring balanced shock absorption. The scooter is equipped with pneumatic tyres and lighting. Furthermore, the tiller is infinitely adjustable, and the armrests can be flipped up and adjusted in width. The Excel Xena EVO is also designed to be dismantled for transport.

Tip

Read this manual in full to ensure optimal use of your mobility scooter.

3. Safety instructions

The Excel Xena EVO mobility scooter is extremely safe and stable, provided that the user follows the instructions outlined in this manual. However, misuse or incorrect operation of the scooter may lead to dangerous situations. Van Os Medical B.V. accepts no responsibility for injury or damage resulting from failure to follow the warnings and recommendations in this manual.

3.1 General safety instructions

Protect your investment by having your scooter checked regularly. If any part of your scooter is not working properly, it may cause a hazardous situation. It is essential to keep your mobility scooter in good condition to ensure safety. Periodic inspections, correct adjustments, and timely replacement of damaged or worn parts will significantly improve the lifespan of the product.

An authorised Excel dealer, using only original Excel spare parts, will help guarantee long-term durability through proper and preventive maintenance.

3.2. Warnings for safe use

To prevent damage to your property or the scooter, or injury to the user, the following warnings must be observed. Van Os Medical B.V. is not liable in any form. By reading these warnings, you accept the terms and conditions as outlined in this user manual. Any unauthorised modifications to the scooter will immediately void the warranty.

- Never drive backwards down a slope.
- Do not drive on uneven surfaces such as sand, gravel or rough terrain, or on wet or slippery surfaces where traction is poor.
- Never drive diagonally up or down a slope.
- Never carry passengers on your scooter. It is designed for single-person use only. Carrying someone on your lap is also not permitted during use.
- Never use your scooter if it is not functioning properly or if any parts are missing or damaged.
- Do not drive up slopes or curbs that exceed the scooter's climbing capability. The maximum safe gradient is 10 degrees.
- Never exceed the stated maximum user weight, which includes bags and any other load on the scooter.
- Never attach anything to the wheels of the scooter. This could cause tipping or damage that compromises safety.
- Never use your scooter on (escalating) stairs, escalators or conveyor belts. See section 6.3.3 for more information.
- Do not use the scooter while under the influence of alcohol, drugs, or medication that may impair your ability to operate it.
- NEVER remain seated in the scooter while it is being transported.
- The scooter can be folded and is suitable for transport in a car boot. Ensure it is securely fastened during transport to avoid movement.
- Transporting the scooter in any vehicle is done at your own risk.

Warning

Making structural modifications to the Excel Xena EVO will void the warranty. Product liability is also forfeited, and all associated risks are transferred to the end user.

- Always check upon delivery whether all components of the scooter are present and undamaged. If anything is missing or damaged, contact your dealer immediately.
- First check that the brakes are functioning properly before setting off.
- Ensure that the anti-tip wheels are installed when driving onto or off of raised surfaces.
- Sudden direction changes may cause imbalance. Steer smoothly, maintain a steady direction while driving, and reduce speed when taking (sharp) turns.
- Always bring the scooter to a complete stop before switching from forward to reverse.
- Ensure the scooter is on a flat surface before getting in or out.
- Always remove the key from the ignition before exiting the scooter.
- Before driving, ensure the seat is facing forward.
- Always use a lift to move the scooter between different levels.
- If the scooter is in a lift or on a scooter lift, it must always be switched off.
- Transporting the scooter in vehicles is entirely at your own risk. Excel Mobility does not accept liability for this.
- Making modifications not described in this manual will void the warranty entirely.

4. Specifications

	Xena EVO	
	3 wheel	4 wheel
Total length (cm)	104	
Total width (cm)	54	
Total height (cm)	88 (excl. mirrors)	
Total weight (kg)	49	52
Max. user weight (kg)	135	
Seat width (cm)	45	
Seat depth (cm)	39	
Seat height (cm)	41	
Ground clearance (cm)	6	
Max. speed (km/h)	6,5	
Range (km)*	17	
Battery	512 Wh	
Charger	72 Wh	

*The driving range is based on manufacturer data measured under ideal conditions. This means: a test weight of 75 kg, 0% battery degradation, 100% battery capacity usage, and a temperature of 25 °C in a controlled indoor environment. Several factors can influence the actual driving range, including user weight, weather conditions, outdoor temperature, speed, tyre pressure, road conditions, battery condition, and the incline of the terrain.

5. Components of your Excel Xena EVO mobility scooter

This chapter describes the various components of your mobility scooter. Several parts are adjustable or removable. If you wish to install, remove, adjust, or reposition any of the components, always follow the instructions provided in this user manual.

5.1 The dashboard

Your mobility scooter is equipped with a control panel. The control panel of the Excel Xena EVO is clearly readable and easy to operate. Photo 2 shows the dashboard, with the components explained below:

Photo 2

Photo 3

1. On/off switch (ignition)

Located at the top of the control panel. Insert the key and turn it 90 degrees clockwise to turn the scooter on. Turn it anticlockwise to switch it off.

Once switched on, nine LED lights will illuminate: three green, three amber, and three red. When all LEDs are lit, the battery is fully charged. As the battery discharges during use, the lights will go out in descending order. If only a single red LED is lit, the battery must be charged. See section 5.6.1 for charging instructions.

2. Horn

Press this button to activate the horn, which emits a beeping sound.

3. Speed control dial

Turn the dial to the right to increase speed or to the left to decrease it.

4. Lighting

The Excel Xena EVO features a front light. Switch it on using the button beneath the dashboard (see photo 4). Use this light to improve visibility in low-light or poor weather conditions.

5. Throttle levers (wigwag system)

To drive forward, pull the right lever towards you (see photo 5).

To reverse, pull the left lever towards you (see photo 6).

The Excel Xena EVO features an automatic braking system: once the lever is released, the scooter will brake automatically.

Photo 4

Photo 5

Photo 6

Photo 7

6. USB port

The dashboard features a USB port (see photo 3), which can be used to connect a navigation-enabled phone, for example.

7. Charging port

Located on the side of the dashboard. See section 5.6.1 for instructions on how to charge the scooter using this port.

5.2 The tiller

The tiller of the Excel Xena EVO can be angle-adjusted to suit your preferred driving position. Follow these steps:

- See photo 8 for the upright position of the tiller.
- On the lower right-hand side of the tiller, there is a black adjustment knob (photo 9). With one hand, loosen the knob while gently pulling the tiller towards you with the other hand (photo 10).
- Once the desired angle is set, retighten the knob.
- To return the tiller to the upright position (recommended for getting out), loosen the knob again and push the tiller away from you.

Photo 8

Photo 9

Photo 10

Photo 11

5.3 The armrests

The Excel Xena EVO features flip-up, width-adjustable, and removable armrests.

To flip them up, simply lift the armrest upwards (photo 12). This makes getting in and out easier (see section 6).

Photo 12

Photo 13

Width adjustment

- Under the seat, there is a black knob (photo 14). To adjust the width, loosen this knob (photo 15) and gently pull the armrest outward (photo 16).
- Once the desired width is set, tighten the knob again.

To completely remove the armrests, fully extend them outward and then pull them off the frame (photo 17).

Photo 14

Photo 15

Photo 16

Photo 17

5.4 The seat

The Excel Xena EVO features a 360-degree swivel seat with a foldable backrest. This makes getting in and out of the scooter easier.

Under the right side of the seat, there is a lever (photo 18). Pulling the lever upwards unlocks the seat, allowing it to rotate. Releasing the lever locks the seat in position.

Photo 18

The upholstery on the seat and backrest is removable:

- Detach the hook-and-loop fastening (Velcro) from the backrest (photo 19).
- Slide the backrest upholstery off (photo 20).
- Repeat the same process for the seat: loosen the Velcro strap at the rear (photo 21) and slide the cover off (photo 22).

Photo 19

Photo 20

Photo 21

Photo 22

5.5 The freewheel

The Excel Xena EVO has two drive modes:

- Drive mode ('D', photo 23)
- Freewheel mode ('N', photo 24)

The mode is indicated on the sticker next to the yellow lever. The scooter must always be in drive mode ('D') when in use.

Freewheel mode allows the scooter to be pushed manually while switched off. In this mode, the motor functions as a generator.

- To activate freewheel mode, pull the yellow lever backwards to the 'N' position.
- To return to drive mode, push the lever forward to 'D' to reactivate the electromagnetic brakes. The scooter can now drive again.

Photo 23

Photo 24

Warning

Before using the Excel Xena EVO, always make sure the freewheel lever is set to drive mode.

5.6 The battery

The Excel Xena EVO is equipped with an LFP (Lithium Iron Phosphate) battery.

Benefits of LFP batteries:

- Higher efficiency for more effective energy use
- Longer lifespan with more charge/discharge cycles
- Higher energy density relative to weight
- Lower maintenance requirements
- Built-in safety features to prevent overcharging, deep discharge, and overheating

5.6.1 Charging the battery

There are two ways to charge the battery. You must follow the steps below to charge the batteries correctly and to extend their lifespan.

Charging via the battery pack

- Can be done at any location because the battery pack is removable (see 5.6.2).
- Always use the standard charger supplied with the scooter.
- Plug the charger into the port on the front of the battery pack (photo 25), then into the socket.
- The LED on the charger will turn amber while charging.
- Once fully charged, the LED will turn green. Disconnect the charger from the mains.
- If the scooter is not in use, recharge the battery once every two months.
- Store the battery in a frost-free room at a minimum of 6 °C.
- LFP batteries can be topped up at any time without shortening their lifespan.

Photo 25

Charging via the dashboard port

You can also charge the batteries via the charging port located on the side of the dashboard. To ensure the batteries are charged correctly and to prolong their lifespan, follow these steps carefully:

- Use only the standard charger supplied with your scooter.
- Plug the charger into the charging port on the side of the dashboard (see photo 26).
- Then insert the other end of the charger into a wall socket.
- An amber LED light on the charger will illuminate to indicate that the batteries are charging.
- Once the batteries are fully charged, a green LED light will appear on the charger. At this point, you must disconnect the charger from the wall socket. You can recognise a fully charged battery by the green LED light illuminated on the charger.

Photo 26

Warning

Only use the charger supplied by Excel Mobility. Using any other charger will void the warranty and may damage the battery.

Warning

Avoid charging the battery in extremely hot or cold temperatures. Use only the original charger. Unauthorised chargers may damage the battery and shorten its life, and will void the warranty.

Tip

Avoid exposure to extreme temperatures, moisture, or direct sunlight. Optimal range is achieved at 25 °C. Below 10 °C, battery performance can drop by up to 30%.

5.6.2 Removing the battery

- The battery is removable (photo 27).
- Grab the handle on the top of the battery (photo 28) and lift it upwards (photo 29) to detach it from the frame.

Photo 27

Photo 28

Photo 29

5.6.3 The fuse

- A thermal fuse is located on the front of the battery (photo 30).
- In the event of a short circuit, the fuse cuts off power to the battery.
- A button on the battery box will pop out. Press it back in to restore power.
- If the button immediately pops out again, the battery is short-circuited—contact your authorised Excel Mobility dealer.

Photo 30

5.7 The anti-tip wheels

The Excel Xena EVO is equipped with rear-mounted anti-tip wheels (photo 31).

These increase safety and help prevent the scooter from tipping backwards—especially when navigating kerbs or slopes.

Photo 31

5.8 The wheels

The Excel Xena EVO is available as a 3-wheel or 4-wheel model. Check the wheels regularly for wear.

5.8.1 Front wheels

- The front wheels are 8" pneumatic tyres (photo 32).
- Check tyre pressure regularly (see 8.1).
- Front wheels are essential for steering. If steering becomes difficult or the wheels vibrate, they may need adjustment.

5.8.2 Rear wheels

- The rear wheels are also 8" pneumatic tyres (photo 33).
- Check these regularly for tyre pressure and wear.

Photo 32

Photo 33

6. User instructions

It is important to familiarise yourself with safe and proper usage techniques. To do so, you must understand all the components described in this user manual. Always consult a professional for support in developing safe usage methods that match your physical abilities.

6.1 Balance

Everyday actions such as getting in and out of the scooter, reaching, and bending will affect the scooter's balance. This has to do with weight distribution and changes to the centre of gravity. To prevent falls, tipping, or other hazardous situations, please follow the guidance below.

6.1.1 Getting in and out

Getting in and out of the mobility scooter is a frequent activity. The level of difficulty may vary depending on your physical condition. This manual assumes that you can enter and exit the scooter independently. If this is no longer the case, we recommend discussing your options with an occupational therapist.

Getting in

- Make sure the scooter is on a flat surface before getting in.
- Always set the scooter to drive mode ('D'). In this mode, the electromagnetic brake is automatically engaged (see section 5.5).
- Ensure the scooter is switched off. The LED lights on the dashboard should be off. This prevents accidental operation of the throttle levers.
- Position the scooter properly for safe access.

You have two options to get in:

- Swivel the seat using the lever on the right side of the seat (turn 90 degrees).
- Flip up the armrest and step in from the side.

Once seated, lower the armrest again and place your feet centrally on the footplate.

If a seatbelt is present, fasten it before driving.

Adjust the tiller angle to a comfortable driving position (see section 5.2), insert the key, and turn it clockwise. The LED lights will come on. You are now ready to drive.

Getting out

Follow the same steps in reverse.

Make sure the scooter is switched off, and return the tiller to the upright position. You can then swivel the seat or flip up the armrest to exit. Try to dismount as close to your destination as possible.

6.2 Bending and/or reaching

You may sometimes need to reach for something while seated. This could be in front of, beside, or behind you. To avoid tipping the scooter or losing your balance, follow these instructions:

6.2.1 Reaching forward

If you want to reach for something in front of your scooter, follow the steps below.

- Drive as close to the object as possible.
- Ensure all wheels are pointing forward.
- Switch off the scooter to engage the brake.
- You can now lean forward safely.

6.2.2 Reaching sideways

If you want to reach for something to the side of your scooter, follow the steps below.

- Drive close to the object, using the rear wheels as a reference point.
- Make sure the wheels are straight and the scooter is switched off.
- Lean to the side with care, keeping your body close to the seat when lifting an object.

6.2.3 Reaching behind

If you want to reach for something behind your scooter, follow the steps below.

- Position the scooter close to the object, preferably against a stable surface such as a cupboard.
- Switch off the scooter to engage the brake.
- Lean back carefully — not too far — to avoid tipping.

6.3 Navigating kerbs and slopes

At times you may need to ascend or descend a kerb or slope. This section explains how to do so safely.

6.3.1 Going up a kerb or slope

- Always approach an incline forwards.
- Drive up at moderate speed — not too slow — to maintain momentum.
- If you stop halfway, restart gently to prevent rollback.

Anti-tip wheels provide added stability when ascending.

Do not hang bags or extra weight on the scooter, as this alters the centre of gravity.

6.3.2 Going down a kerb or slope

- Always descend forwards and in a straight line.
- Reduce speed for maximum control.
- Use the throttle gently, and keep the scooter in its lowest speed setting.

Tip

Always inspect the slope or kerb before driving over it so you know what to expect.

Warning

If the scooter speeds up more than you're comfortable with, release the throttle slowly. The scooter will automatically brake and reduce speed.

Warning

Never attempt to drive on or off slopes covered in snow, ice, wet grass, cut grass, leaves, or other surfaces that may reduce grip.

6.3.3 Moving surfaces

The Excel Xena EVO is not suitable for moving surfaces such as escalators or conveyor belts. As a powered vehicle, it should only be used on stable, solid surfaces. If available, always use a lift.

6.4 Steering

Steering is one of the most important aspects of driving a mobility scooter. Below you will find instructions for safe steering in various situations.

- Always place both hands on the tiller grips when steering.
- Keep your arms slightly bent rather than fully extended — this improves comfort and control.
- When turning, always reduce your speed to avoid dangerous situations or tipping.

6.4.1 Steering in different situations

The Excel Xena EVO is designed to enhance your mobility, and you may encounter a variety of environments. Below are a few common situations and how to handle them safely:

Steering in tight spaces

- Bring the scooter to a complete stop.
- Set the speed control to the lowest setting.
- Turn the tiller in the desired direction.
- Gently engage the throttle lever.

Steering in reverse

- Set the speed to low.
- Use your left hand to pull the reverse throttle lever.
- Turn the tiller left to go left, or right to go right.
- If your scooter is fitted with mirrors, use them when reversing.

Challenging surfaces

- The scooter performs best on flat surfaces.
- It is possible to ride on uneven terrain such as grass, gravel, or sand, or even smooth icy surfaces — but this is not recommended.
- If you do so, reduce your speed to match the conditions.
- Avoid coming to a stop on such surfaces by keeping the scooter in the lowest speed setting.

6.5 EMI / RFI

With the rapid development of communication technologies, our surroundings are increasingly filled with electromagnetic (radio) waves — emitted by TV transmitters, mobile phones, radios, laptops, microwave ovens, and pacemakers. These waves are invisible and become stronger the closer you are to the source.

When electromagnetic waves disrupt or cause uncontrolled behaviour in electronic devices, this is known as Electromagnetic Interference (EMI) or Radio Frequency Interference (RFI).

All electrical vehicles, including mobility scooters, are vulnerable to EMI/RFI, which may cause unexpected or abnormal movements. If this happens — for example, if the scooter behaves erratically or the brakes do not work — switch the scooter off immediately.

The FDA has specified that all mobility scooters must withstand EMI/RFI up to a minimum level. This level is measured in volts per metre (V/m). The minimum resistance is 20 V/m. This product has been tested and complies with that standard.

6.5.1 EMI / RFI recommendations

- Do not use personal electronic communication devices such as mobile phones while operating the scooter.
- Be especially cautious around transmission masts and radio/TV towers.
- Switch off your mobility scooter when stationary for extended periods.

7 Transport and handling of the mobility scooter

The Excel Xena EVO mobility scooter has been designed with easy transportation in mind. The scooter can be dismantled for compact storage and transport in a car.

7.1 Dismantling the scooter for transport

To dismantle the scooter quickly and easily, follow the steps below:

- Remove the key from the ignition.
- Remove the front basket (photo 34) by lifting it upwards.
- Remove the seat by pulling the lever on the right-hand side (photo 35), gripping both sides of the seat, and lifting it upwards.
- You may remove the armrests to reduce the weight (see section 5.3).
- Remove the battery pack from the frame (see section 5.6.2 and photo 36).
- Fold down the tiller by loosening the black adjustment knob located underneath on the right side (photo 37). Then fold the tiller down and retighten the knob.
- Lock the tiller in place by pressing the black locking knob (photo 38) and turning it a quarter turn clockwise (photo 39).
- Separate the front and rear frames by lifting the black release handle and gently pulling the rear frame away from the front frame (photo 40).
- The scooter is now fully dismantled (photo 41) and ready for transport.

Photo 34

Photo 35

Photo 36

Photo 37

Photo 38

Photo 39

Photo 40

Photo 41

7.2 Reassembling the scooter for use

To reassemble the scooter, follow these steps:

- Reconnect the rear frame to the front frame by lifting the black release handle and placing the rear section into the designated yellow slots. The frame will click into place.
- Unlock the tiller by pressing the black locking knob (photo 42) and turning it a quarter turn anticlockwise (photo 43).
- Raise the tiller into the upright position by loosening the adjustment knob (photo 44), lifting the tiller, and tightening the knob again.
- Reattach the front basket to the tiller.
- Place the battery pack back into the frame.
- Reinstall the seat by gripping both sides and lowering it into the seat post tube (photo 45).
- If you removed the armrests, reattach them (see section 5.3).
- Insert the key into the ignition, and the scooter is ready for use.

Photo 42

Photo 43

Photo 44

Photo 45

8 Maintenance

Your Excel Xena EVO mobility scooter requires regular maintenance. This is essential for long-term durability and optimal user comfort. A poorly maintained scooter may experience technical issues, run less smoothly, and fall outside of warranty coverage. Preventive maintenance is therefore highly recommended.

We advise having your scooter inspected at least once a year by an authorised Excel Mobility dealer. The dealer will use only genuine Excel replacement parts when necessary. With annual servicing, your scooter will continue to perform reliably for many years.

You can also do a lot yourself to keep your scooter in top condition. By performing regular checks and carrying out minor maintenance when needed, you will prolong its lifespan and improve your everyday experience. The following sections describe what you can do yourself.

8.1 Tyres

- Check the tyres regularly for wear. If they are significantly worn, they must be replaced.
- Pneumatic tyres must be checked weekly for proper pressure.
- The correct pressure is indicated on the tyre itself.
- Use a standard pressure gauge to check.
- If the pressure is incorrect, the scooter may lean to one side, which increases wear, impairs steering, and requires more effort to drive.
- Clean the tyres with car shampoo and a damp cloth. Do not use solvents, as they may damage or soften the material.

8.2 Bodywork

Proper maintenance of the bodywork is important, as mechanical parts like the motor often outlast an unmaintained frame. The goal is to prevent rust for as long as possible.

- Wash the scooter regularly (once per week recommended) to prevent dirt build-up that can damage the paint.
- Clean the bodywork with car shampoo and a damp cloth.
- Do not spray the scooter with a pressure washer.
- Dry the scooter with a soft cloth.
- Use cold water and mild soap to remove dirt and grease.

8.3 Upholstery

The seat upholstery is also an important safety element. Torn or worn upholstery can lead to unsafe situations.

- Check the upholstery regularly for holes, tears, or signs of wear.
- If you notice any damage, contact your dealer. Only an authorised dealer may perform repairs to the upholstery.

8.4 Cleaning your mobility scooter

Regular cleaning is important. The best way to clean your scooter is as follows:

- Clean the upholstery, chassis, and plastic parts with mild soap and water.
- Do not use abrasive cleaners, as they may damage the paint.
- Never use a pressure washer or steam cleaner.
- Polish the painted parts of the chassis regularly using car wax.
- Avoid solvent-based or abrasive wax, harmful chemicals, or silicone sprays.
- Always dry the scooter thoroughly after cleaning, or after exposure to rain.

8.5 Periodic checks

Preventive maintenance is important, and many tasks can be done by yourself or with help from a caregiver. We recommend the following checklist:

Every week

- Check the tyre pressure.

Every month

- Inspect the upholstery for damage or wear.
- Check that the armrests fold up smoothly. Lubricate the joints with silicone spray if needed.
- Ensure the seat rotates easily. If not, lubricate the mechanism with silicone spray.

Every quarter

- Check the tread on the tyres for signs of wear.

9. Troubleshooting

If you experience issues while using your Excel mobility scooter, refer to the table below for possible causes and solutions. If your issue is not listed, we recommend contacting an authorised Excel Mobility dealer.

Number of LED flashes	Solution
1x	The batteries need charging or there is a poor connection to the battery. Check all connections and plug in the charger.
2x	Battery voltage is too low. Charge the batteries and inspect both the batteries and wiring.
3x	Battery voltage is too high. If going downhill, reduce speed to limit regenerative charging.
4x	Current limit timeout / controller overheated. The motor has drawn maximum current for too long. The scooter may stop. Turn off the controller, wait a few minutes, then switch back on. Check the motor and wiring.
5x	Brake fault. Check the parking brake and its connections and wiring. Make sure all switches are in the correct position.
6x	Stop function activated, charging lock-out, or OONAPU condition. Release the stop condition. Disconnect the charger and ensure the throttle lever is in neutral before turning the controller back on.
7x	Speed potentiometer fault. Check the potentiometer and related connections and wiring.
8x	Brake connection fault. Check connections between the brake and motor. Ensure all control system connections are correctly attached.
9x	Other internal fault. Contact an authorised dealer.

10. Warranty

Your purchase invoice serves as your proof of warranty.

10.1 Warranty coverage

Your mobility scooter is covered by a manufacturer's warranty. In some cases, your supplier may offer an additional warranty, but this section covers only the factory warranty provided by Excel Mobility.

The manufacturer's warranty applies only to private users. It is not valid for commercial or institutional use (e.g. rentals or care facilities). The warranty is limited to material defects and hidden faults. Excel Mobility offers the following warranty periods:

- Complete product: 2 years
- EVO battery: 5 years
- EVO charger: 2 years

10.2 Warranty terms and conditions

The warranty period starts on the date of purchase. To claim a warranty repair, you must return the scooter to your dealer and provide the original purchase invoice.

Repairs and replacements are carried out by authorised Excel Mobility dealers.

If the scooter is used in a different area from the original dealer, the repair may be handled by another authorised service partner.

If a component covered under warranty is found to be defective due to a manufacturing or material fault — and the scooter is still owned by the original purchaser — the part will be repaired or replaced free of charge by the authorised dealer.

Note: Labour costs for fitting or replacing parts are not covered under warranty.

Repaired or replaced parts are subject to the same warranty conditions as the original product. Wear-and-tear items are not normally covered, unless failure is due to a factory fault.

Excel Mobility is not responsible for repairs or replacements caused by:

- Failure to maintain the scooter according to the manufacturer's instructions.
- Use of non-original parts or components.
- Damage caused by negligence, accidents, or misuse.
- Unauthorised modifications or repairs not carried out by an authorised service agent.
- Missing or removed factory frame labels (see section 1).

The mobility scooter described and shown in this user manual may differ in detail from your actual model. However, all instructions are relevant regardless of these differences. We reserve the right to make changes to the weight, dimensions, or other technical specifications described in this manual without prior notice. All illustrations, dimensions, and specifications shown in this user manual are approximate and may not fully correspond to the exact features of your model.

For warranty service, please contact the authorised dealer from whom you purchased your mobility scooter. If you are not satisfied with the assistance provided by your dealer, you may contact Excel Mobility directly. The importer's contact details can be found on the back of this manual.

Warning

Excel Mobility cannot be held liable for any incidental or consequential damages. This manual has been written with care, but is not exhaustive. If you intend to use the scooter in a way not described in this manual, consult an authorised Excel Mobility dealer first.

Warning

Warranty is valid only within the specified period. Any structural modifications to the scooter will void the warranty entirely. Contact Excel Mobility for complete warranty and delivery conditions, or for a list of authorised dealers.

Tip

Maintenance is not included under warranty. Your dealer may advise a different service interval.

10.2.1 Warranty terms for LiFePO4 batteries

1. a. Battery warranty period: The warranty period for LiFePO4 batteries is 5 years or 2000 charge cycles, whichever occurs first.
b. Charger warranty period: The warranty period for the lithium battery charger is 2 years from the date of purchase.
2. Warranty claims: During the warranty period, Excel Mobility will repair or replace the LiFePO4 battery and/or charger free of charge if a defect occurs due to a manufacturing or material fault, or if battery degradation exceeds the limits stated in the table below.
3. Conditions for warranty validity: To make a warranty claim, the owner must be able to demonstrate that the battery and/or charger has been used, charged, and discharged according to the manufacturer's recommendations, and in accordance with the instructions provided in the user manual.
4. Warranty exclusions: The warranty does not cover any defects or damage caused by improper use or maintenance, use of an unauthorised charger, unauthorised modifications, misuse, accidents, negligence, or external factors beyond Excel Mobility's control.
5. Warranty limitations: The warranty is valid only for the original owner of the scooter and is non-transferable.
6. Battery replacement: If the battery and/or charger needs to be replaced under warranty, Excel Mobility will provide a new or equivalent battery or charger, at their discretion.
7. Liability: Excel Mobility is not liable for any indirect, incidental, or consequential damage arising from the use — or inability to use — the battery and/or charger.
8. Claim procedure: To submit a warranty claim for a LiFePO4 battery and/or charger, contact your authorised Excel Mobility dealer. They will provide the necessary instructions and support to handle the claim.

For more information about the warranty terms for LiFePO4 batteries, please contact Excel Mobility.

The following degradation percentages apply to your battery. Degradation is calculated based on the battery's original net capacity:

Year	Annual degradation	Cumulative degradation	Charge cycles
1	4%	4%	300
2	3%	7%	600
3	2%	9%	900
4	2%	11%	1200
5	2%	13%	1500
6	1%	14%	1800
7	1%	15%	2100

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